


Graphic designer/ p. lotfolahi Photo By M.cheraghchi

SINA

NEW PRODUCT



Salimi ind Town, Azarshahr road,
Tabriz-Iran P.O.Box:53751-131
www.sinabed.com
Tel: +98(41)34 32 92 42
Fax: +98(41)34 32 92 43

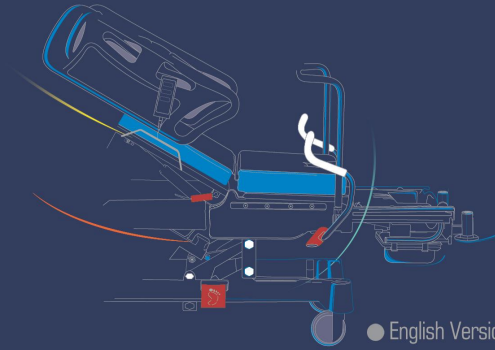

● English Version 2018 / CASTIA

Castia Series

Sina newly designed delivery beds that ease the trauma of birth by providing comfort and support throughout the stages of delivery. A selection of accessories increases the comfort of both the mother and nursing staff.


This kind of bed helps to keep the state of a natural childbirth as healthy as possible. This multi-functional bed due to the special characteristics, provides comfort feeling for pregnant women before, during and after delivery by catching three position. The mother is admitted to the same bed throughout the delivery and early days after childbirth. In this way, many of the problems and discomforts caused by the transfer of mother to stretcher or labor is removed.



Castia Series
D 3101 poc

Height and trend . / Rev . trend and backrest are adjusted by means of 3 actuators and 3 - Function handset . CPR handle is available at both sides of the CASTIA for flattening bed rapidly . 4- segment mattress is washable and disinfected . Footboard and leg - section mattress are removed in lithotomic position. Those are placed behind the headboard . At this position leg - section is placed under the seat .

Castia Series
D3101 pose

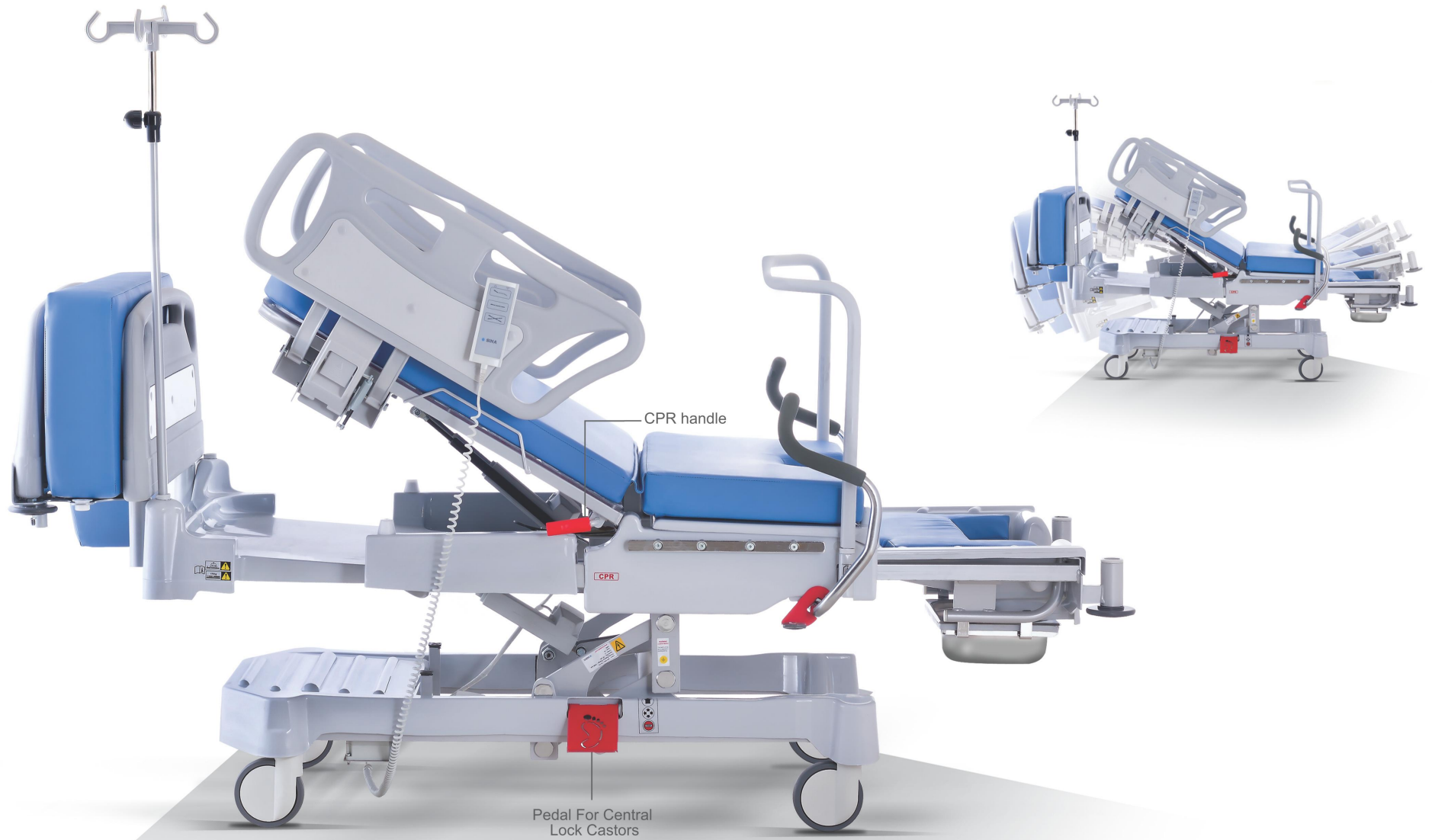

Castia Series



Plastic sideboards are completely safe and userfriendly. Squatting bar and hand grip device are available when needed. In case of need, leg - section mattress is removed, then underlay mattress appears. All parts of CASTIA are washable and disinfectable.

TECHNICAL INFORMATION

- Backrest, Hi- Low and Trend/ Rev. Trend adjustments are electrically
- Two segmented plastic sideboard.
- Three- function positioning control accessible to the patient.
- With ABS cover for covering all parts.
- Holder for squatting position.
- Leg section can be placed under seat for lithotomic positions
- pull out and removable bowl (5 litre)
- with hand grip in both sides
- Detachable head/ footboards
- Washable and disinfectable
- CPR handle on both sides for flattening bed rapidly
- Backup battery
- Four section detachable mattress.
- Four castors with central/directional brakes
- Four rubber bumpers at each corner.
- Main power supply: 220V, 50 Hz, 24 DC
- Max. Power input: 500 watt
- Electrical current protection: Class1
- Electrical shock position: Type B
- Duty cycle: %10
- Dust proof and water-resistant: IP66
- IV stand holder at each corner.
- Traction holder.
- Monkey Bar holder.

DIMENSIONS AND BASIC INFORMATION

- Over all: 2340 X 1000 mm
- Matters platform: 2000 X 865mm.
- Back rest angle: 0 – 80
- Hi-low: 595 – 785 mm
- Trend/Rev.Trend:+ 15 /- 15
- Castors: 150 mm diameter central/ directional brakes
- Max. load: 250kg.

ACCESSORIES

- IV stand
- Monkey Bar
- Traction Bar