

SINA

Atria

TECHNICAL CATALOGUE



سینا متعدد

Atria





Atria is a new mechanical hospital bed with a unique design. The strength platform has a distinct design. This puzzle shaped platform includes plastic segments which is aligned on a chassis. SINA's recent innovation in adjusting the Backrest and Knee raise is PUSH - PUSH Cranks which is user friendly. These special cranks with freewheeling system warn the user when are in Max. Min angle. Special ergonomic designed Head/ Foot board is another features. Plastofold side rails are completely safe and user friendly.




● PUSH - PUSH crank for kneeraise and backrest adjustment



● Traction holder



● Puzzle shaped platform

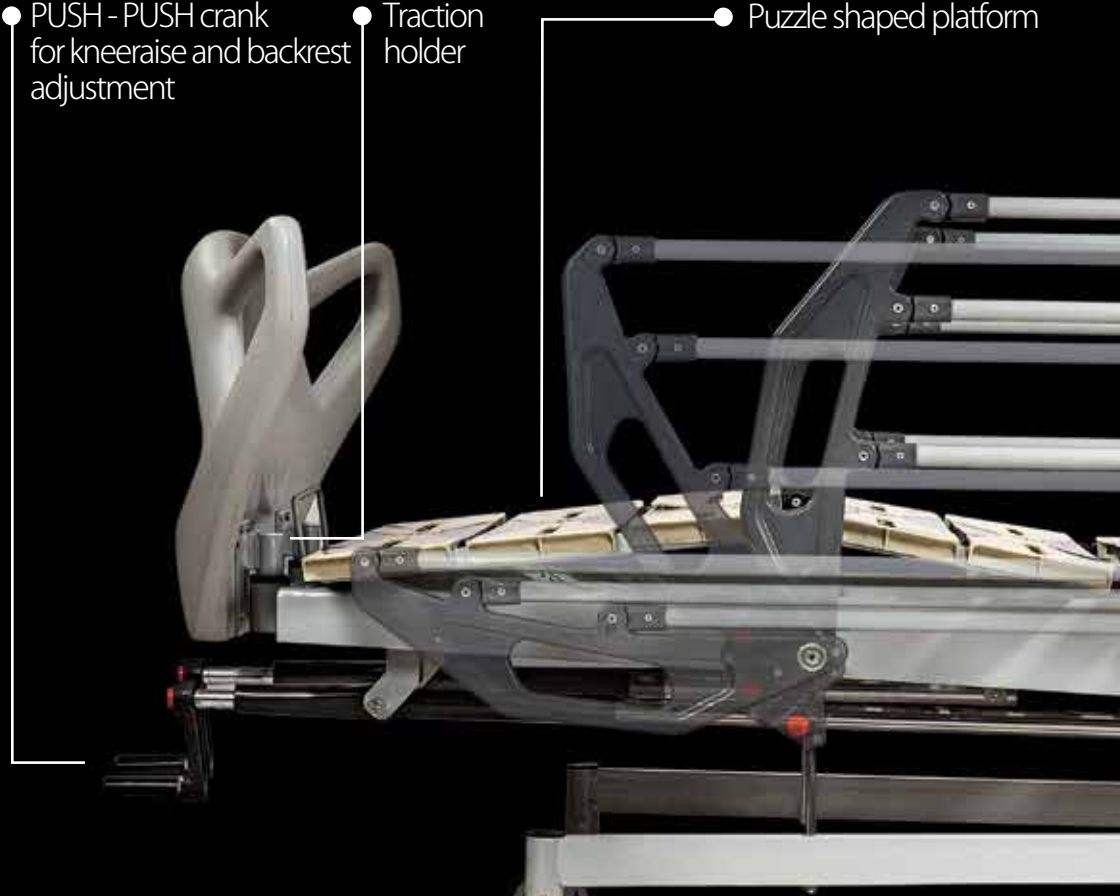


● Independent break castors



● Plastofold sideboard



● Plastic head / foot board with metal connection

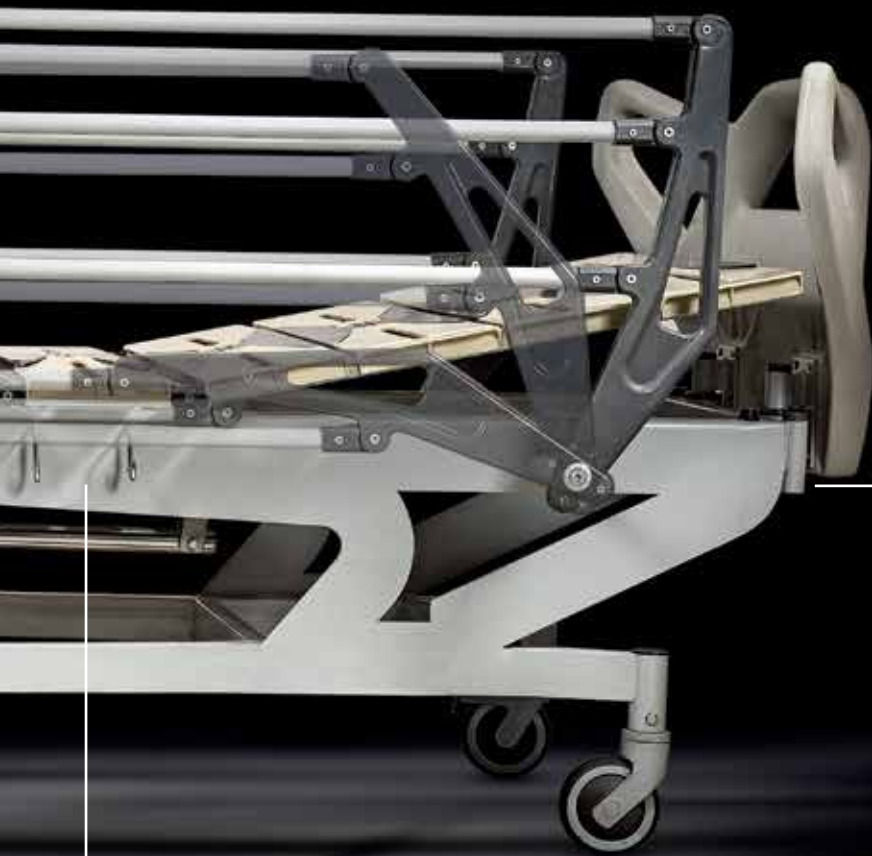

● IV Stand holder



● Urine bag holder



TECHNICAL INFORMATION

- Mattress platform is four section, three cranks with puzzle shaped system and strengthened chassis
- Strong hinges made of polyamide
- Plastofold sideboards
- Foldable PUSH - PUSH cranks with freewheeling system
- Washable and disinfectable
- Air flow capability under the mattress
- Four casters with 2 independent brake system
- Plastic Head/foot board with metal connection and suitable ergonomic style for easy handling
- With IV Stand in four corners of the bed
- Urine bag holders are available at both side of the bed
- Traction holder is available at foot board







DIMENSION and BASIC INFORMATION

- Overall 229 x 102 x 59 (cm)
- Mattress platform 200 x 90 (cm)
- Backrest angle 0 - 85
- Knee raise angle 0 - 35
- Casters 125mm
- Weight 90 kg
- Max. load 250kg

Standards

- 60601-1 IEC
- 60601-2-38 IEC
- ISO 9001
- ISO 13485



SINA

20m st. of 33, North 3rd
30m st., 1st 45m st. Salimi
Ind. town, Tabriz - Azarshahr
road, Iran - Tabriz
www.sinabed.com
info@sinabed.com
Tel: +98(41) 34 32 92 42
Fax: +98(41) 34 32 92 43