

TECHNICAL CATALOGUE



سینا متعدد

Stia





Stra is a new product in stretchers group. Stra's remarkable design is based on being user friendly and high safety for users.

The unique features of Stra are:

- Appropriate hydraulic lifts
- Linear bumpers along the stretcher
- Special transfer handle and smooth castors



Height and Trend/Rev. Trend is adjustable by the set of 3 pedals. These pedals are available under the cover of stretcher, with a set on each side. This feature gives the nurse freedom to use his/her hands for patient control and at the same time changing the position with foot quickly without wasting any time. Using levers on side rails instead of a button has made it so convenient.

Stra is appropriate for using in emergency ward, Operating room, recovery ward, etc.







TECHNICAL INFORMATION

- Mattress platform is in two section with one crack
- Detachable Mattress
- Hi-Low and Trend/Rev. Trend is adjustable by hydraulic lifts
- Backrest is adjustable by gas spring
- Leverside Sideboards
- Washable and disinfect able
- 200 cm diameter castors with central/directional brakes.
- With bumpers in four corners of the bed and the linear bumpers in both side of the stretcher
- Oxygen cylinder holder on the cover of the bed
- I.V stand holder is available at four corners
- With a holder for patient file



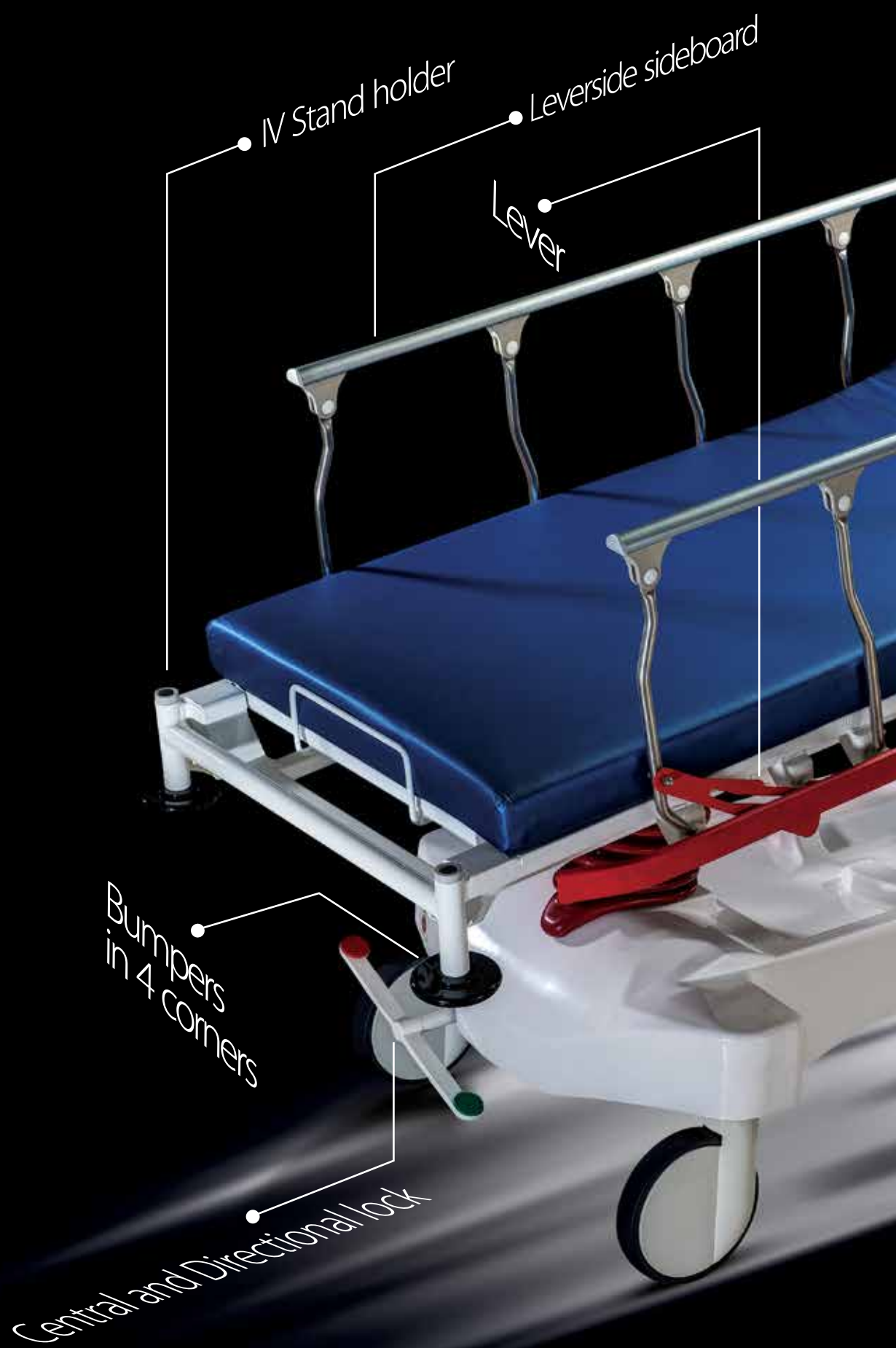

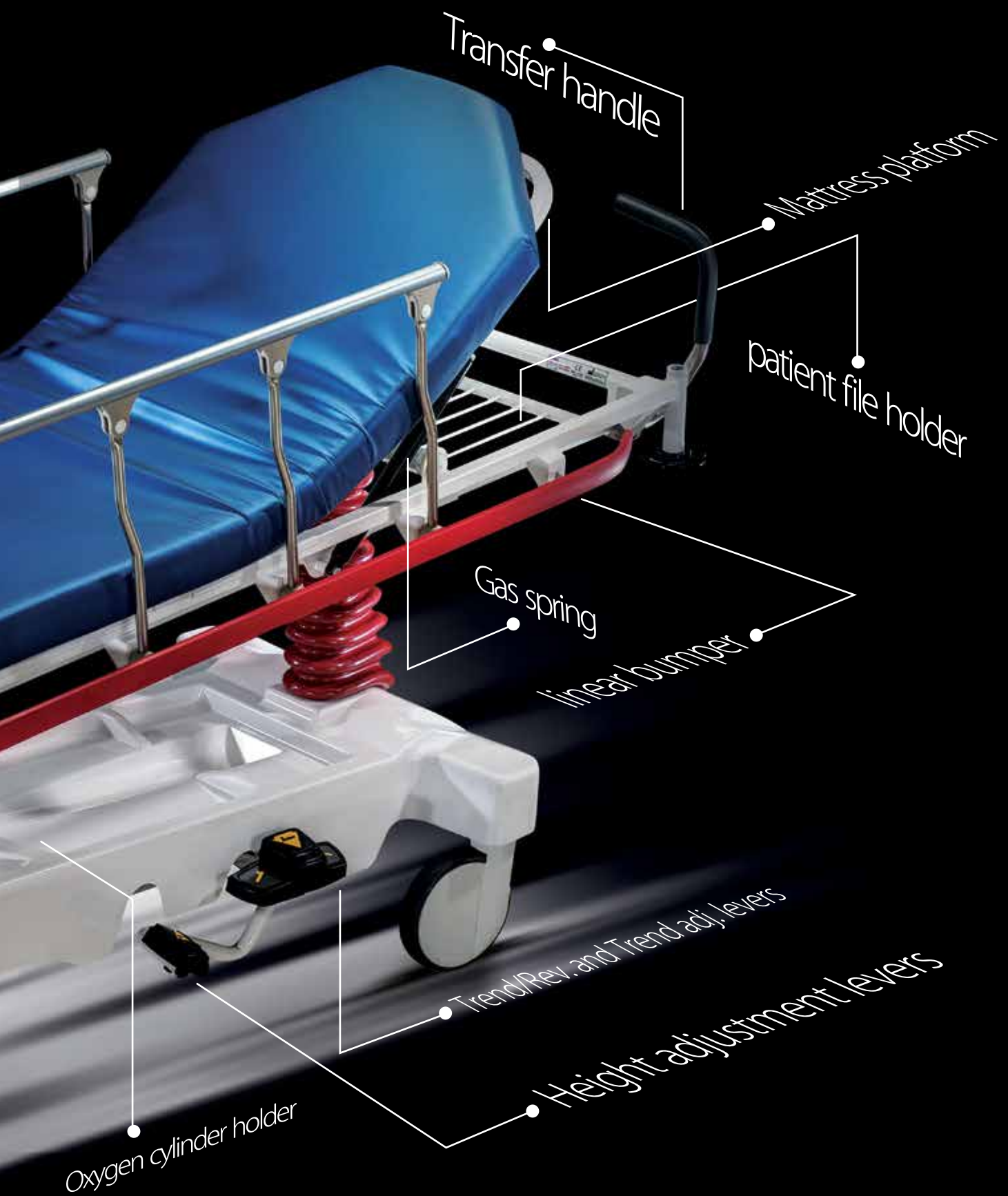






DIMENSIONS and BASIC INFORMATION

- Overall 2114 x 754 x 870 (mm)
- Mattress platform 1865 x 635 (mm)
- Hi-Low 550 - 870 (mm)
- Backrest angle 0° - 90°
- Trend / Rev. Trend angle +15°/-15°
- Castors 200 mm
- Weight 94 kg
- Max. load 250 kg

Standards

- IEC- 60601-1
- IEC- 60601-2-52
- ISO 9001
- ISO 13485

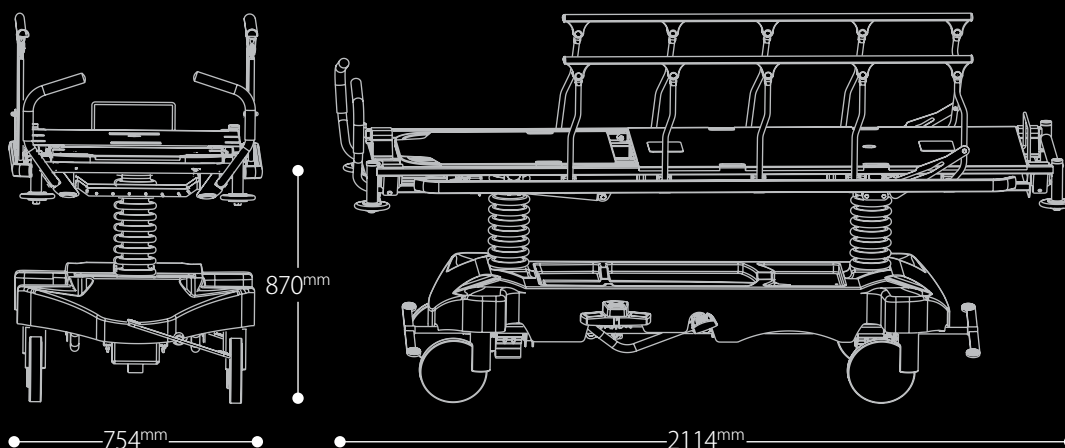

SINA

20m st. of 33, North 3rd
30m st., 1st 45m st. Salimi
Ind. town, Tabriz - Azarshahr
road, Iran - Tabriz
www.sinabed.com
info@sinabed.com
Tel: +98(41) 34 32 92 42
Fax: +98(41) 34 32 92 43

