


Hepatopathy

BS330 Plasma Adsorption Cartridge (For hyperbilirubinemia and hyperbileacidemia)

BS330, filling with anion-exchange resin, can specifically **adsorb bilirubin** and **bile acids**. It works effectively to relieve the symptoms of **hyperbilirubinemia** and **hyperbileacidemia**.

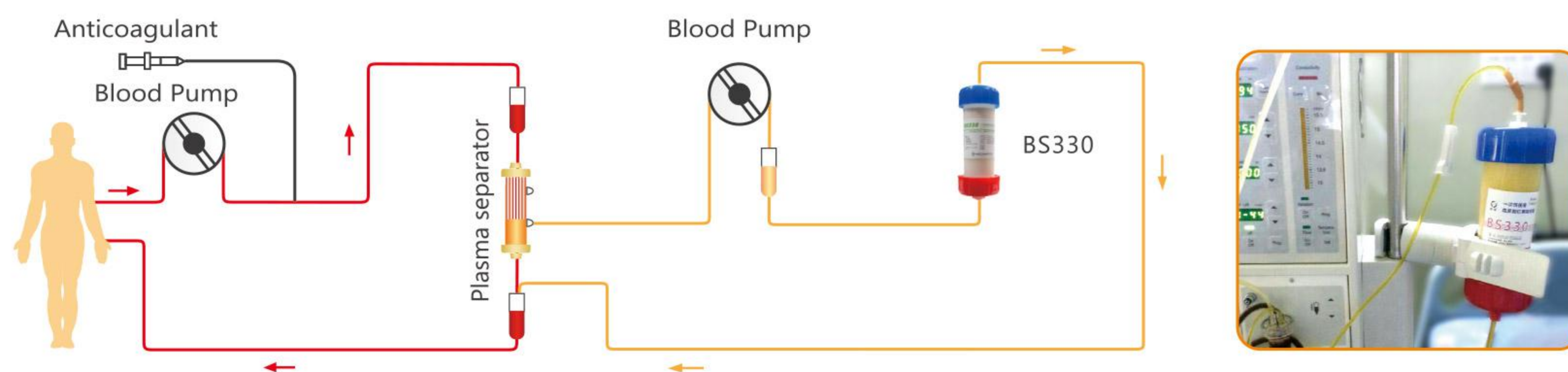


Laboratory Test Results


Adsorption test result of plasma in vitro, the efficacy of adsorption for 6 hours:

adsorption duration (h)	circulating plasma volume	the average decreased total bilirubin
2	3600	59.68%
4	7200	67.16%
6	10800	71.06%

Clinic Trial of BS330

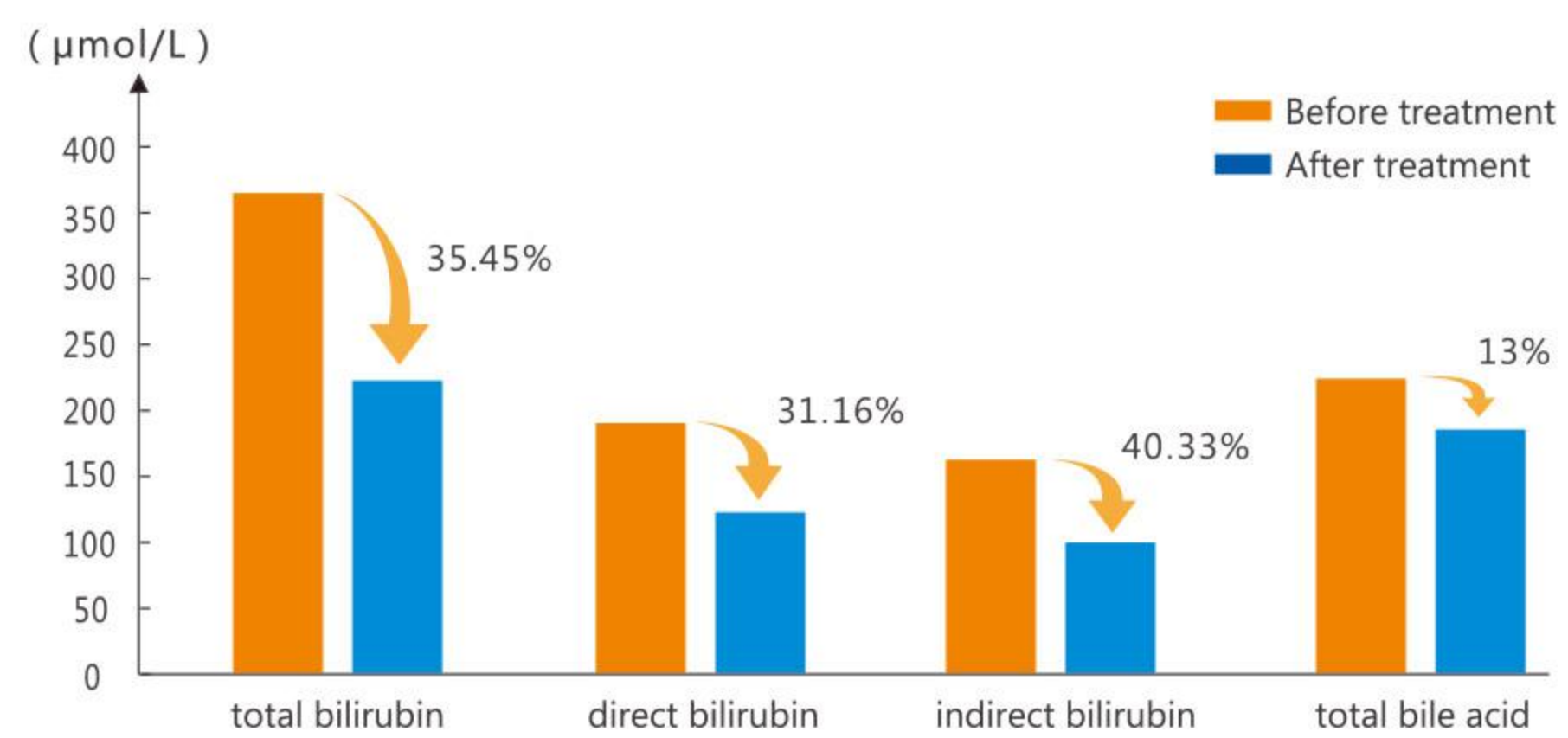

BS330 + HA330- II (DPMAS) Comprehensive Solution for Hepatopathy

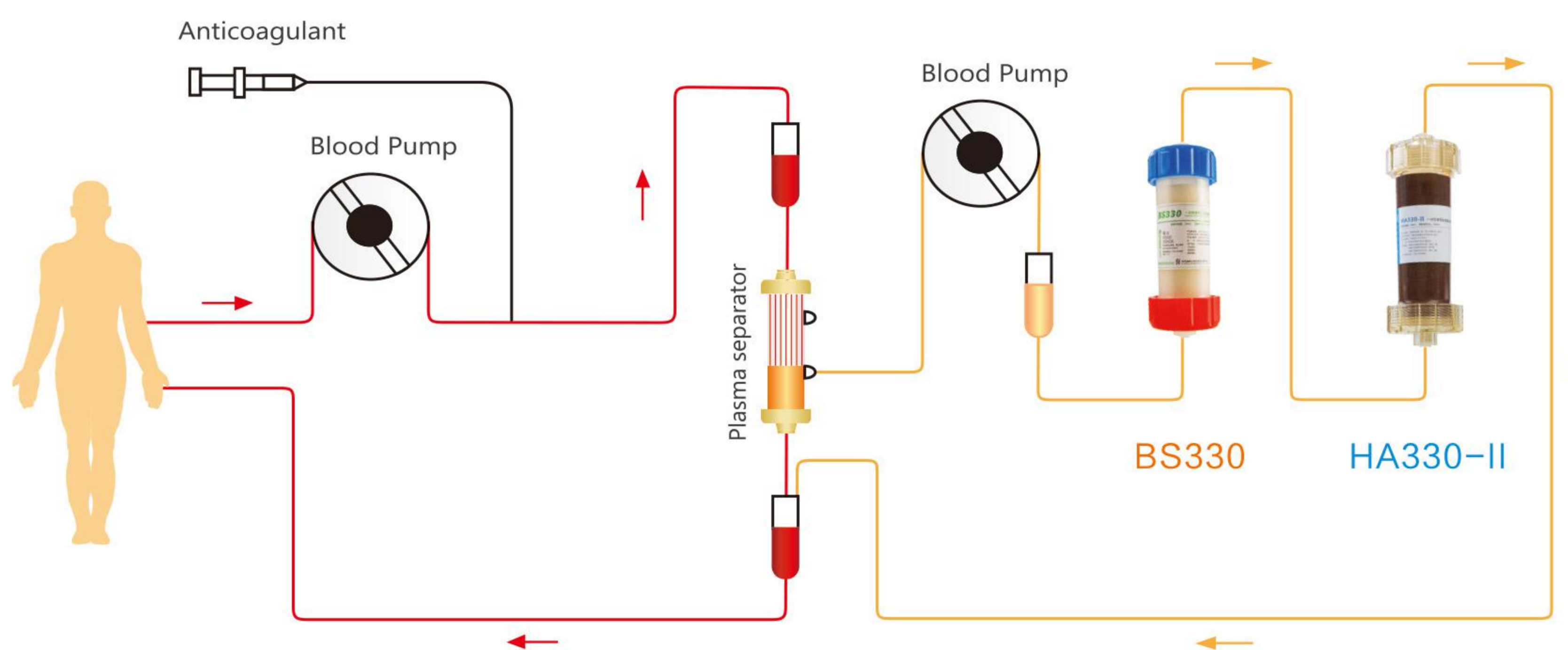


HA330- II Disposable Hemoperfusion Cartridge

(For liver disease)

HA330- II extensively adsorbs **middle molecular toxins** induced by **liver disorder**, such as inflammatory mediators, ammonia, phenol, mercaptan, etc.

Recommend for drug-induced hepatitis, hepatic failure & complications such as hepatic encephalopathy, etc.

Flexible Operation

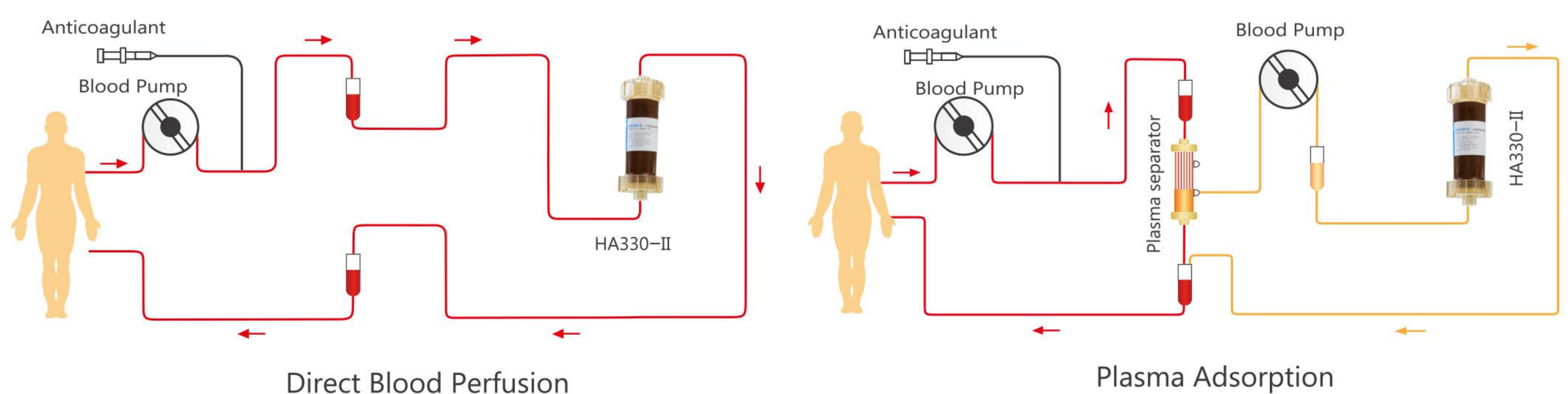

The Double Plasma Molecular Adsorption System (DPMAS) specifically adsorbs **bilirubin**, effectively clears **inflammatory mediators**, eases inflammation and immune responses as well as considerably relieves the symptoms of hepatopathy and eventually improves the **long-term prognosis** of patients.

